

Transform Connectivity

OneLab Secured Deployment Solution



OneLab




design & execute

Seamless installation, instant functionality, and effortless integration —combined with your lab's scientific expertise—create the perfect recipe for breakthroughs.



OneLab™ Software puts scientists in control with intuitive protocol design, powerful automation, and a user-friendly graphical interface that provides access to an extensive online protocol library. OneLab Secured Deployment Solution elevates this experience with a secure, turnkey Windows® Operating System-based workstation that integrates seamlessly with your network, ensures cybersecurity compliance, and simplifies IT management for hassle-free deployment.

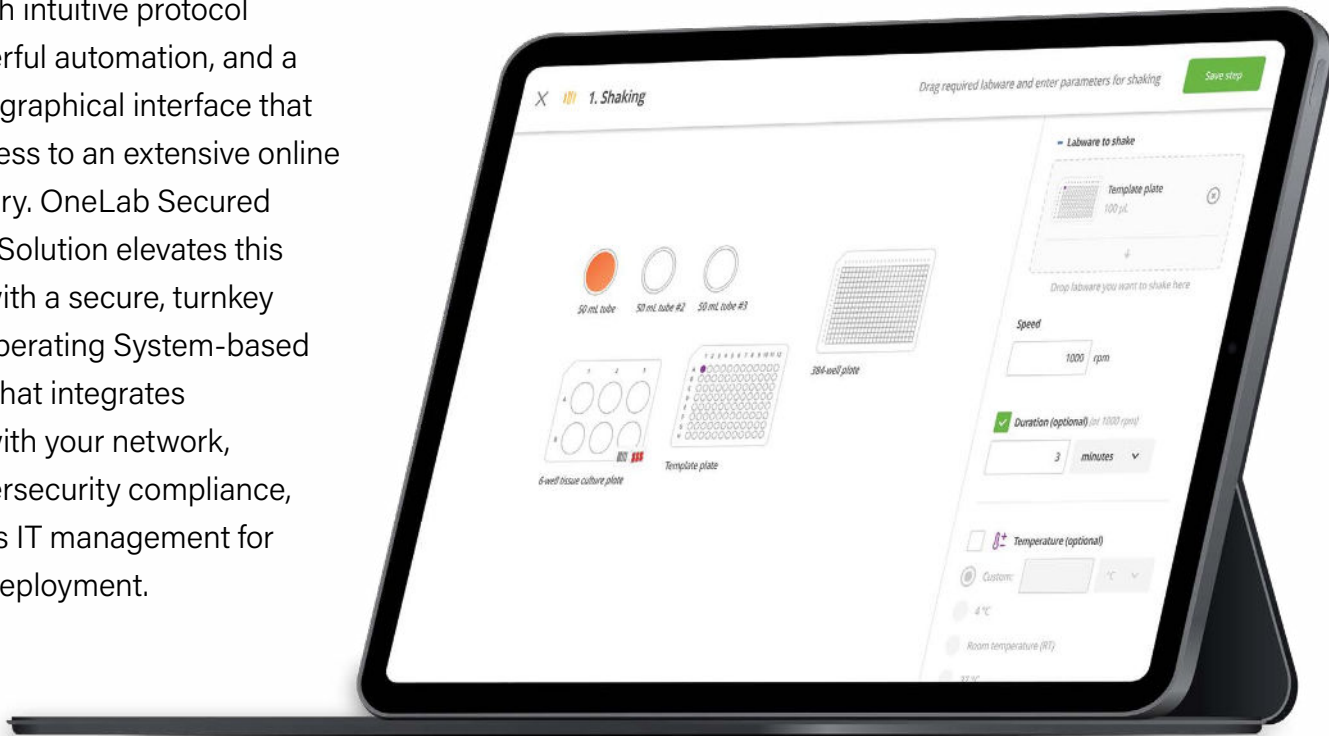

Box to Bench in a Blink

Connect, install, and start automating in minutes—without complex IT setup.



Enterprise-Ready Performance

Enable smooth operation, effortless updates and maintenance, and secure backup options with validated hardware and software.



Uncompromising Cybersecurity

Meet strict enterprise security standards for safe operation in regulated environments.



Compliance You Can Count On

Access compliant-ready features that support 21 CFR Part 11 and GxP needs, ensuring confidence and readiness for regulated environments.




Future-Ready Integrations

Grows with your lab and aligns with your IT infrastructure to deliver secure, seamless integrations for long-term scalability.


Waters Validated Platform

Waters carefully selects and qualifies hardware partners to guarantee consistent quality, global availability, and an optimized footprint for laboratory environments. Each OneLab Secured Deployment Solution workstation is rigorously tested to ensure seamless compatibility with all versions of OneLab Software, delivering reliability and compliance-readiness without compromise.

OneLab Secured Deployment Solution Device Overview

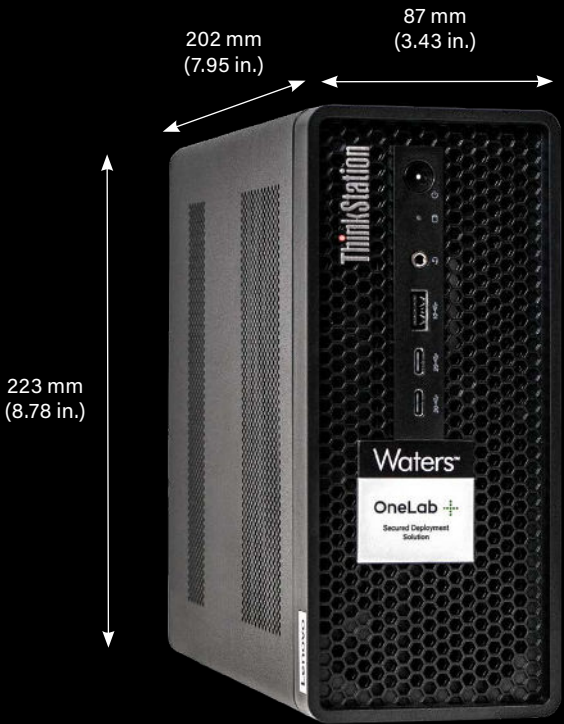

Based on Microsoft® Windows Enterprise LTSC (Long-Term Servicing Channel)

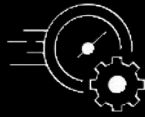

Compact chassis and premium performance




One Ethernet port for corporate network connection



Up to six Ethernet ports for connected devices and or user computer

- Isolated from the corporate network
- An Ethernet port is required for each connected device
- This local network can be expanded by using an additional Ethernet switch

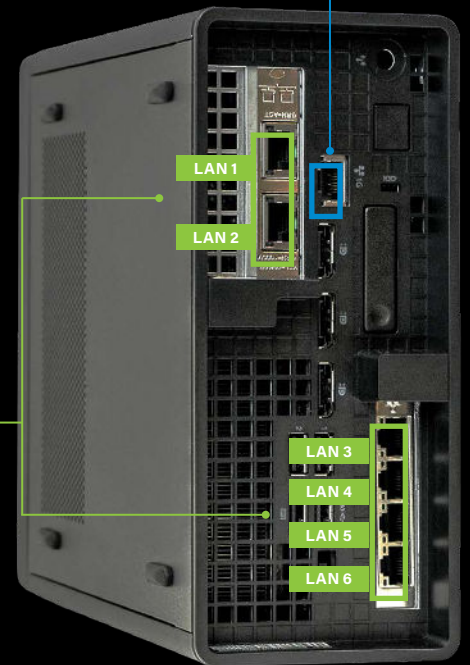


Connected for Smarter Science

- Unlock true interoperability with Empower™ Software, MassLynx™ MS Software, and waters_connect™ Software solutions
- Automate high-quality data capture from OneLab Software workflows
- Streamline your processes so scientists can focus on results, not routine tasks



Your Digital Lab Console

Connect a screen to the **OneLab Secured Deployment Solution** workstation and use it as a console for your system. A keyboard and mouse are included with every workstation.

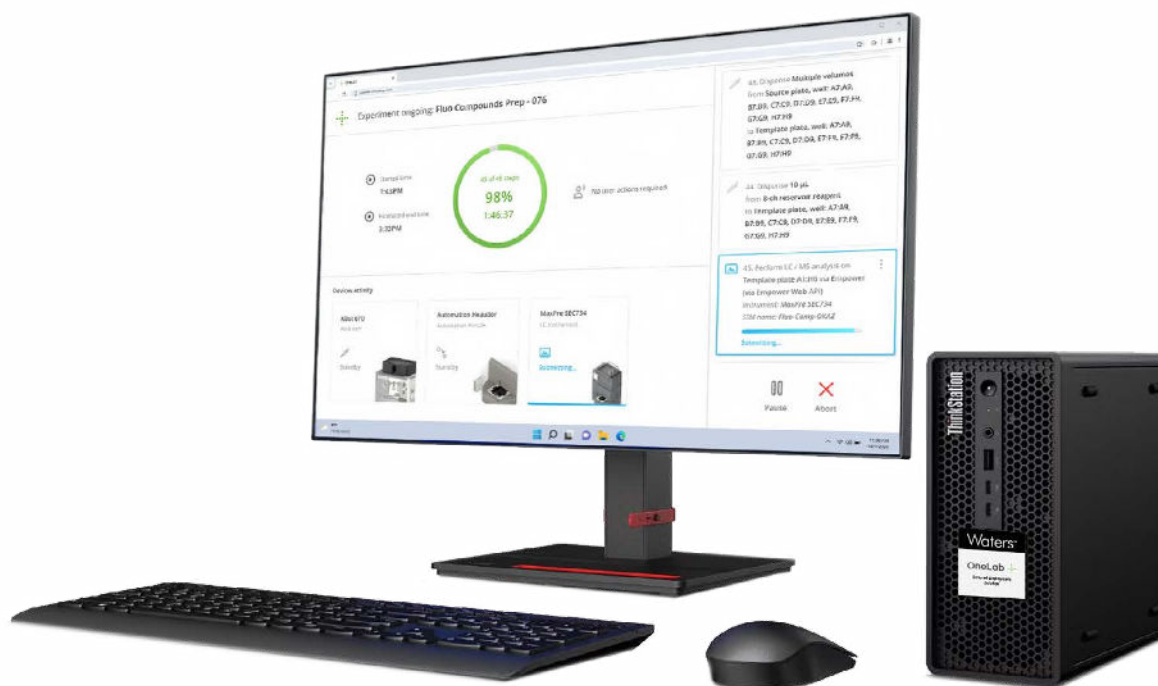

Ordering Information

All-in-One OneLab Secured Deployment Solution SDS Bundle

Complete, ready-to-deploy package, including:

- Preconfigured Windows 11 LTSC computer
- Keyboard, mouse & DisplayPort-HDMI adapter
- OneLab SDS Installer
- OneLab Enterprise 10 users License (1-year)

Bundle Options

- 1-Year — P/N 176006219
- 3-Year — P/N 176006220
- 5-Year — P/N 176006221

Licenses for Expansion or Renewal

10 users OneLab Enterprise licenses for OneLab SDS:

- 1-Year — P/N 667007906DD
- 3-Year — P/N 667007907DD
- 5-Year — P/N 667007908DD

Computer monitor option

- Lenovo 24" ThinkVision Flat Panel: P/N 668000273

Technical Specifications

Operating System	Windows 11 Enterprise LTSC
Processor	Intel® Core™ Ultra 5 235 (E-cores up to 4.40 GHz, P-cores up to 5.00 GHz)
Memory	32 GB DDR5-5600 MT/s (SODIMM)
Storage	512 GB SSD M.2 2280 PCIe Gen4 Performance TLC Opal
Graphics	Integrated Graphics, DisplayPort connector
Ethernet	I350-T4 (4x Ports PCIe 1G LP)
Second Ethernet	BCM5720 Dual-Port PCIe 1G Adapter
Keyboard	USB, Traditional, Black - English (EU)
Mouse	USB Calliope Mouse (Black with Red Wheel)
External Power Supply Input Characteristics	Input: 100–240 VAC, 4.5 A, 50–60 Hz Output: 20 VDC, 16.5 A, 330 W
Warranty	4Y Premier Support Plus (Upgrade from 3Y Onsite)
Country of Manufacture	Hungary
Dimensions (mm) W x D x H	202 mm x 87 mm x 223 mm
Weight	1.8 Kg

waters.com/[labautomation](https://waters.com/labautomation)

For your local sales
office, please visit
waters.com/contact



Waters Corporation
34 Maple Street
Milford, MA 01757 U.S.A.
T: 1 508 478 2000
F: 1 508 872 1990
waters.com

Waters™

Waters, OneLab, Andrew+, Pipette+, MassLynx, Empower, and waters_connect are
trademarks of Waters Technologies Corporation.
All other trademarks are the property of their respective owners.

©2026 Waters Corporation. January 25-14850 720009204EN REV. A